

# matrx-pb

POSTURE BACK

The Matrx-PB Posture Backs are contoured to provide midline support and positioning. The wide variety of sizes and depths will provide exceptional comfort and support for a wide range of user needs.



## BACK:

- Contoured shell (standard/deep) provides exceptional midline support and positioning
- Foam wraps around shell for maximum comfort and protection
- Durable, lightweight ABS back shell includes built-in plate for headrest attachment
- Moisture resistant, breathable cover
- Includes removable lumbar positioning pad

## HARDWARE:

- "Easy Set" hardware installed and adjusted in minutes – adjustments can be made while user is seated in the wheelchair
- Quick release for easy portability
- Infinite adjustment of depth (4"), angle (40°) and rotation
- Fits 3/4", 7/8" and 1" diameter back canes
- Adjusts to fit chairs 1" narrower to 1" wider



### STANDARD POSTURE BACK & DEEP POSTURE BACK SIZES

width	14"	16"	17"	18"	20"
height					
12"	●	●	●	●	●
16"	●	●	●	●	●
20"		●	●	●	●

### HEAVY DUTY POSTURE BACK & DEEP HEAVY DUTY POSTURE BACK SIZES

width	20"	22"	24"
height			
16"	■	■	■
20"	■	■	■

### EXTRA WIDE HEAVY DUTY POSTURE BACK & DEEP EXTRA WIDE HEAVY DUTY POSTURE BACK SIZES

width	24"	25"	26"	27"	28"	29"	30"
height							
16"	□	□	□	□	□	□	□
20"	□	□	□	□	□	□	□



Kid\*ab\*ra fabric optional

● Weight capacity 300 lbs.

■ Weight capacity 500 lbs. □ Weight capacity 500 lbs.

## PRODUCT OPTIONS

### FOAM-IN-PLACE KIT

- Individually molded insert produced in minutes with user seated in back
- Optional high stretch cover available

### EXTRA LONG MOUNTING PINS




- Extended mounting pins available for Matrx-PB
- Allows the back to fit chairs up to 2" wider

### VISCO FOAM UPGRADE

- Laminate of 1/2" high resilient foam and 5/8" Visco Foam
- Available on Matrx PB and Elite backs

### CANE REPLACEMENT

- Replaces existing back canes which are less than 1" diameter for use with any Matrx back

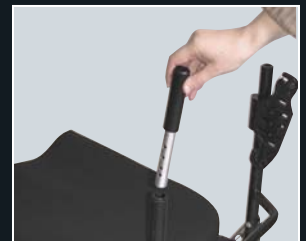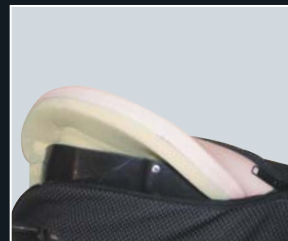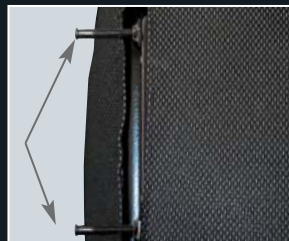





**Standard "Easy Set" PB Mounting Hardware (PBMH)**



- Included with standard MaTRx PB Posture/Elite backs
- Infinite adjustment of depth (4"), recline (30°) and even rotation
- Fits 7/8" or 1" diameter back canes (inserts available for 3/4" back)

**Heavy Duty "Easy Set" Mounting Hardware (HDMH)**



- Standard with MaTRx PB Heavy Duty/Posture backs
- Recommended for users over 300lbs, users with very high tone or those who are very hard on equipment

**Split Mounting Hardware (SMH)**



- For use with standard MaTRx PB Posture Elite backs - not recommended for Heavy Duty backs
- Provides clearance when mounting on wheelchairs with rigidizer bars midway on the back canes

**Fixed Space Saver Mounting Hardware (FSSH)**



- Narrows mounting hardware by 1" for tighter width applications
- Non-portable - tools required for removal
- Allows rotation, depth and angle adjustment in ESR an PSB interfaces

**Extended Shear Reduction Mounting Hardware**



- Allows Elite backs and Posture backs to be mounted directly to TRx power seating

**Posture Back Cane Replacement**



- Replaces existing back cane extensions which are less than 1" diameter
- For use with any MaTRx back

**Posture Back Cane Clamp Inserts**



- 7/8" inserts included with all MaTRx PB Posture/Elite backs
- 3/4" inserts no charge with any MaTRx back